


Home Oxygen Equipment And Consumables Paediatrics



baywater
H E A L T H C A R E

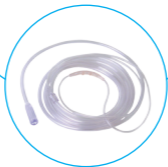

Paediatric Cannula Flexicare

Premature

0-4months, 0.6lpm

Neonatal

0-9months, 0-6lpm

Infant

9months-5years, 0-6lpm

Paediatric

6years-12 years, 0-6lpm



Paediatric Cannula **Intersurgical**

Neonatal


0-4lpm

Infant (Phthalatate free)

0-4lpm

Paediatric

0-4lpm

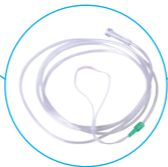

Paediatric Cannula **Salter**

Newborn/Neonatal
Newborn, 0-3lpm

Intermediate
9months-1year, 0-3lpm

Infant
9months-5years, 0-3lpm

High Flow
9months-5years, 0-6lpm



Paediatric Cannula Salter (continued)

Paediatric

6-12years, 0-6lpm

Straight Adult

12years and over, 0-5lpm

Curved Adult

12years and over, 0-5lpm

High Flow Adult

12years and over, 6-8lpm

Recommended Care of Disposables

Nasal Cannula

Wipe tip of cannula daily with a damp cloth. Do not immerse in water.

Dispose of cannula monthly or after swimming/bathing or dirty.

Tubing

Change only if kinked, damaged or discoloured.

It is recommended to dispose of any mask or cannula after an infection.



Venturi Mask Oxygen Concentration

24%: 2lpm **Blue**

28%: 2-4lpm **White**

31%: 6lpm **Orange**

35%: 8lpm **Yellow**

40%: 10lpm **Red**

60%: 15lpm **Green**



Other Masks

Paediatric Tracheostomy
Mask (Venturi only)

Paediatric Eco Mask (soft)
24%-60% (Venturi only)

Paediatric Medium
Concentration Mask
2-8lpm

Paediatric Non
Re-breathe Mask 8-15lpm

Recommended Care of Disposables

Masks

Wash mask in hot soapy water daily and dispose of 3 – 6 monthly, if the mask is visibly dirty or marked.

Caution if using a non rebreathe mask, do not immerse in water. Clean mask with a damp cloth only.

Tubing

Change only if very kinked or discoloured.



For More Information

Clinician Hotline:

01270 218050

Portal:

<https://goo.gl/XscMoN>

Website:

www.baywater.co.uk

